

# EZee Life™

---

## Tilt Overbed Table

Model: CH2006



EZee Life™ Products

34 Futurity Gate, Unit 15  
Vaughan, Ontario L4K 1S6

## Limited Warranty:

Ezee Life™ products carry a one year warranty against manufacturing defects, faulty materials and workmanship only.

Any alterations, misuse, abuse, or accidental damage voids this warranty. Ezee Life™ reserves the right to replace or repair any part that has become defective.

The warranty does not extend to consequential costs resulting from the fault or defect of the product: freight and travel costs, loss of earnings, or other expenses that one may incur.

Warranty is valid for the original purchaser only and the original receipt must be presented before any warranty options are considered. Please also be prepared to provide the date of purchase and serial number.

## Specifications:

Large table section width:	23.5"
Small table section width:	8"
Table depth:	15"
Overall table top size:	32 3/4" wide x 15" deep
Table height range:	29" to 49"
Construction:	Steel base / Laminate table top
Product weight:	19.4 lbs

## Features:

- Swivel casters for easy manoeuvrability
- Locking casters
- Angle adjustable table top
- Beveled edges (top and bottom)
- Easy clean table top
- Large top section adjusts to any angle
- Small top section stays horizontal
- Height adjustable

## Assembly:

The CH2006 overbed table is delivered in three simple pieces and the necessary hardware for assembly.

1. Table top
2. Wheeled base
3. Adjustable table post



Step 1: Align table post with two holes in base ensuring the tab is positioned over the hole in the caster bar. Insert screws and tighten.



Table Post Bottom

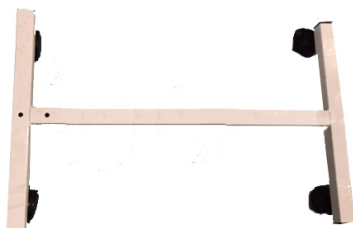


Table Base

Step 2: Align four holes in table post top plate with the four holes in table top bracket. Ensure the base is positioned under the table. Insert four bolts, screw nuts on and tighten.



Table Post Top



Table Top Bracket

## Adjusting Table Height

Support the table top with one hand while loosening the clamp knob (4) on the top of the table post.

Loosen the clamp knob and adjust table top height. Retighten clamp knob ensuring the table top is secured in position.



## Adjusting Table Angle

To adjust the overbed table top angle loosen the clamp lever (5) and pivot the large section of the table top to the desired angle. Retighten the clamp lever ensuring the table angle is secure.

**Note:** The table top has ridges on both sides to prevent items such as books and papers from sliding off when tilted.

**Note:** The small section of the table top does not tilt so there is a level place for drinks or other spillable items.

## Locking Casters

The CH2006 is equipped with four swivel locking casters. The casters can be locked and unlocked by pressing the lock tabs on each caster down or lifting them up using one's toe.

**Note:** Locking and unlocking the casters will be easier and more comfortable when wearing footwear.

### Imported / Distributed by:

Factory Direct Medical  
34 Futurity Gate, Unit # 15  
Concord Ontario L4K 1S6  
Local: 416-739-8393  
Toll Free: 855-235-2400